

STAGE IV®
MILLENNIUM PLUS

STAGE IV®
3000 PLUS

THERA-TURN®
MILLENNIUM PLUS

Alternating Pressure
Low-Air-Loss
Lateral Rotation
Bariatric Mattress
Replacement Systems



Building Distinction and
Excellence into Bariatric
Patient Care

Superior **comfort** and advanced wound care **therapy** now available for the **bariatric** patient.

AVAILABLE
FOR
BARIATRICS




SUPERIOR
ALTERNATING
PRESSURE



ADJUSTABLE
AIR
FLOTATION



LOW AIR LOSS
MICROVENT
TECHNOLOGY

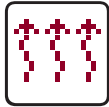

MAXIMUM
SHEAR
RELIEF




The problem of obesity is at epidemic proportions with more than 14 million people in the United States seriously overweight and about 5 million so seriously overweight it affects health and life expectancy. According to the National Institutes of Health obesity is now the second leading cause of preventable death in the United States.

The treatment and care of bariatric patients presents special challenges for providers and caregivers. In particular, bariatric patients with limited or no mobility are at extremely high risk for pressure ulcer development. Products must be specially designed to meet the needs of these patients while providing optimum therapeutic results.

SenTech Medical, a leader in specialty support surfaces, has designed and engineered a full product line of support surfaces to meet special bariatric needs. SenTech's team of engineers, working closely with caregivers, is proud to offer the most advanced line of support surfaces designed specifically for the bariatric patient.

The STAGE IV PLUS MATTRESS

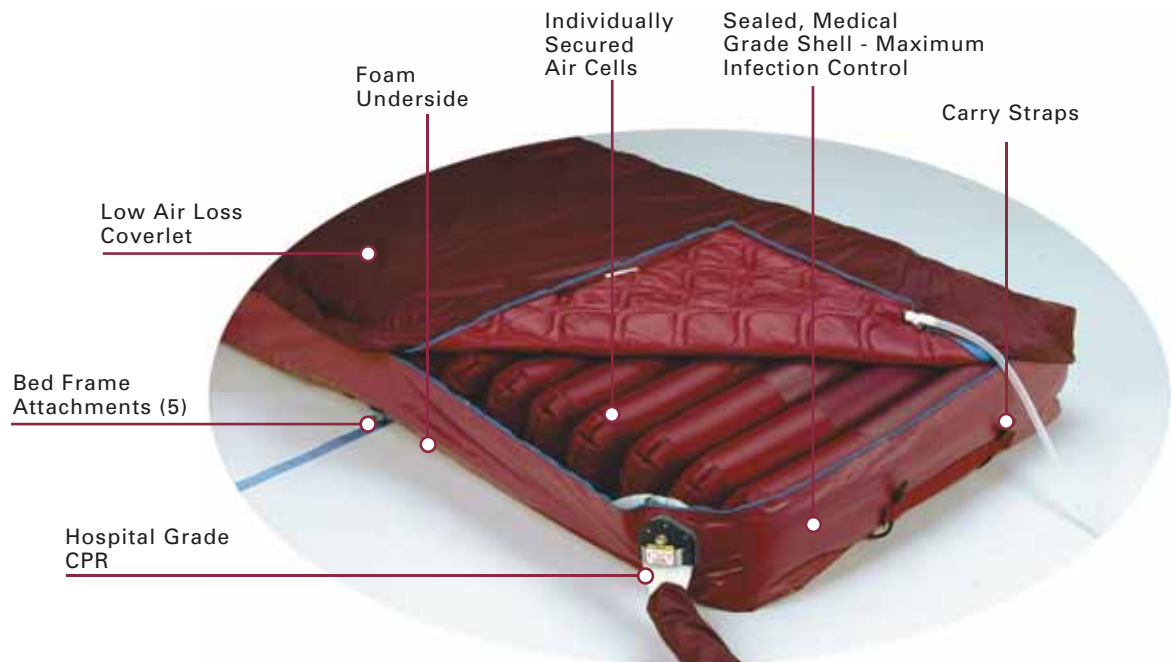

STAGE IV BARIATRICS SYSTEMS

The STAGE IV 3000 PLUS and STAGE IV MILLENNIUM PLUS

Supports patients up to 990 lbs.



Stage IV 3000 PLUS:

- Multi-Zoned System (the head, trunk and foot sections each can be individually customized)
- Sensor Technology (Internal sensors detect and adjust changes in pressure)
- Built-in Memory (stored settings, even after a power loss).
- Visual Monitoring of Therapy (Caregivers receive continuous feedback of therapy readings)
- Patient Height and Weight Settings
- Panel Lock-out feature
- Fowler Boost
- Deflate mode



The Stage IV 3000 PLUS multi-zoned system.



The Stage IV Millennium PLUS as Easy as 1, 2, 3.

- 1 Enter Patient Height
- 2 Enter Patient Weight
- 3 Press Auto Set



The Stage IV Millennium PLUS system provides continuous interactive feedback.

Stage IV Millennium PLUS:

Everything our Stage IV 3000 PLUS Controller has, and...

- Simple to use Interactive Set-up Screen
- Internal Operating Instructions
- Context Specific Help Screen
- Adjustable Max Firm, Fowler Boost Pressures
- Cycle Time Adjustments from 4-9 minutes
- Alarm Locator (displays area to check if alarm is activated)
- Usage Meter (displays hours of usage)
- Filter Display (automatically reminds you to change filter)

Unparalleled, State of the Art, Lateral Rotation

Now Available for the Bariatric Patient.

AVAILABLE FOR BARIATRICS



3 LAYERED AIR CELL PLUS SYSTEM




LOW AIR LOSS MICROVENT TECHNOLOGY



MAXIMUM SHEAR RELIEF




HIGH DEGREE OF TURN

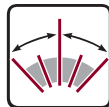

SenTech's new Thera-Turn Millennium Plus is the most sophisticated, yet easy to use lateral rotation system on the market today. And now it's available for the Bariatric Patient.

Incorporating state of the art digital and sensor technology, while supporting patient weights up to 750 pounds, the Thera-Turn Millennium Plus provides superior rotation along with unmatched patient comfort and unparalleled flexibility.

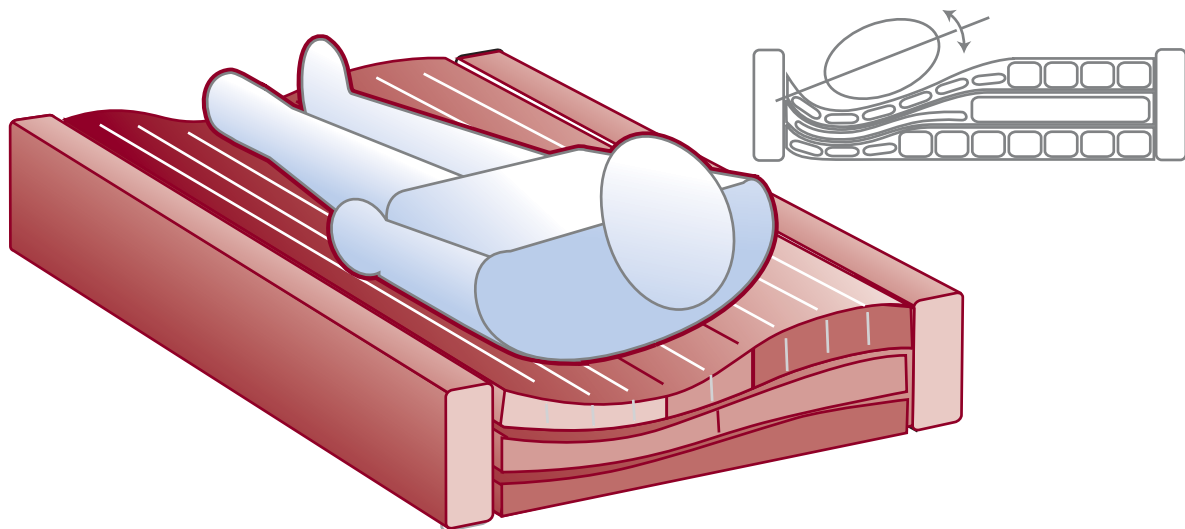

Thera-Turn Millennium Plus Controller



- Easy Interactive Setup Screen
- 4 Pre-Programmed Turning Angles
- Choice of Hold Positions Center, Left or Right
- Choice of Angle-of-Hold
- Comfort Level Setting
- Automatic Filter Display
- Internal Operating Instructions

Specifications Stage IV Millennium Plus & Stage IV 3000 Plus

Electronic Controller:

Dimensions: 14.5" W x 6" D x 12.1"H
 Weight: 21 lbs
 Electrical Parameters:
 Power Requirements- 120V, 60 Cycle AC
 Europe: 240V, 50 Cycle <0.5A
 Maximum Current-<1.0A
 Operating Parameters:
 Weight Range:
 3000 Plus-80-990 lbs.
 Millennium Plus-80-995 lbs.
 Height Range-50"-78"

Mattress Replacement (fully inflated):

Height (Inches) 8
 Width: (Inches) 42,48,54 & 60
 Length: (Inches) 80
 Weight: 24-45 lbs
 Custom Sizes Available

Specifications Thera-Turn Millennium Plus

Electronic Controller:

Dimensions: 14.5" W x 6" D x 12.1"H
 Weight: 21 lbs
 Electrical Parameters:
 Power Requirements- 120V, 60 Cycle AC
 Europe: 240V, 50 Cycle <0.5A
 Maximum Current-<1.0A
 Operating Parameters:
 Weight Range-250 -750 lbs.
 Height Range-50"-75"
 Turn Time Range- 4-90 min. (left/right) 0-90 (center).

Mattress Replacement (fully inflated):

Dimensions: 48"W x 80"D x 14" H
 Weight: 42 lbs

SenTech and Thera-Turn Millennium Plus are Trademarks of SenTech Medical Systems Inc. Stage IV, Millennium Plus and MicroVent Technology are Registered Trademarks of SenTech Medical Systems Inc.

SenTechTM
 Medical Systems, Inc.

5353 N.W. 35th Ave.
 Fort Lauderdale, Florida 33309
 Phone: 1.800.474.4225
 954.739.4972
 Fax: 954.739.2779
www.sentechmedical.com
 REV 40904