

CARROLL CS SERIES™
CS9™ FX600
ADJUSTABLE-WIDTH BED



CS9 FX600 bed featured with Invacare® ThinkSoft® Positioning Device, Kendall Bed Ends and Invacare® Solace® Mattress

The innovative **Carroll CS Series™ CS9™ FX600 Adjustable Width Bed** adjusts from 36" to 39" or 42" to accommodate a variety of residents including those up to 600 lb., offering facilities the utmost in flexibility and providing a total solution for any long-term care facility. Matching the look and feel of the entire CS line, this FX600 is a bariatric bed that blends in with the others, bringing dignity to the resident. The optional trendlenburg option makes this bed a true solution in every wing of the facility.

Features

Expandable Width Capability – Adjusts from 36" to 39" to 42" without additional parts or tools. Tool-less adjustable headboard and 39" fixed footboard keep the dimensions within the entrapment guidelines.

Advanced Electronics – Synchronized motors designed to keep the bed level at all times, providing stability throughout the travel range.

One-Touch Floor Lock – Patented one-touch floor lock secures the bed at any height.



Pure Vertical Rise – Raises straight up and down throughout the travel range, with zero horizontal movement of the casters minimizing potential damage to the wall or flooring.

Exclusive Slat Deck Design – Designed to provide a comfortable sleep surface, allows for improved infection control and reduces wear and tear on the mattress.

Ergonomic Hand Pendant – Waterproof, rugged hand pendant designed to fit comfortably in the resident's hand with dual pendant connections.

Customizable Options –

- Attendant Control Panel
- Trendelenburg / Reverse Trendelenburg / Comfort Chair
- Battery Backup on board or portable
- Hi-Lo Lockout
- 4 Rail options

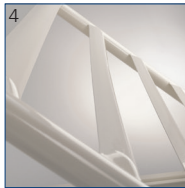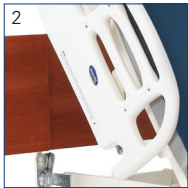

1. Expandable Deck Width
2. Expandable Headboard
3. One-Touch Floor Lock

4. Slat Deck Design
5. 6-Function Ergonomic Hand Pendant

Invacare Corporation
www.invacare-cc.com

USA

One Invacare Way
Elyria, Ohio
44035-4190
(800) 668-2337

Canada

570 Matheson Blvd. E.,
Unit 8
Mississauga, Ontario
L4Z 4G4 Canada
(800) 668-5324

©2016

Invacare Corporation.
All rights reserved.
Trademarks are identified by the symbols™ and ©. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted. Carroll CS Series and CS9 are trademarks of Carroll Healthcare. Specifications are subject to change without notification. Form No. 12-636 Rev. 02/16 150618



Model #	Height Range	Sleep Surface	Weight Capacity*	Regulation Compliance	Limited Warranty
IHCS9FX	9" - 30" H	76" or 80" L 36" - 39" - 42" W	600 lb.	UL 60601-1 IEC 60601-2-38 CAN/CSA C22.2 No. 601.1	10 Years Frame 5 Years Mechanical 3 Years Electronics

*Maximum safe working load.

